

Week 6: Electrocardiography (ECG) II

BIOE 320 Systems Physiology Laboratory

Objectives

1. To record ECG measurements from leads I and III in the following conditions: lying down, sitting up, and breathing deeply while sitting.
2. To use Einthoven's Law to estimate the mean electrical axis and magnitude of the QRS complex.

Background

Electrical Activity of the Heart (Review)

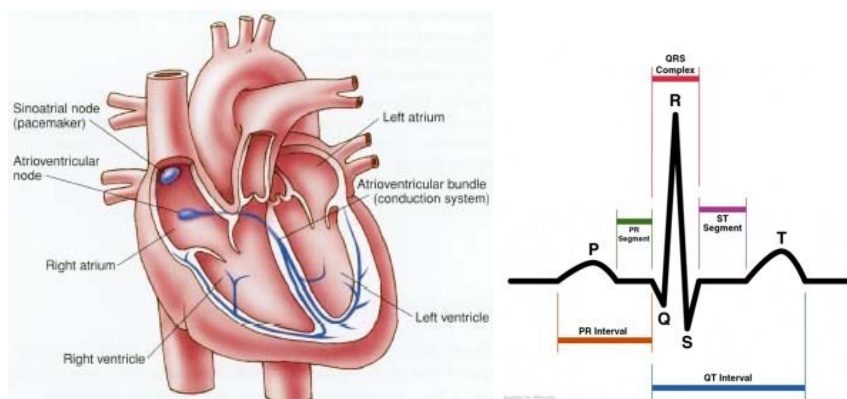


Figure 1: The pacemaker system and the PQRST wave

The heart's electrical system begins at the sinoatrial (SA) node, the primary pacemaker of the heart, where electric potential originates. The polarization spreads throughout the atria, initiating contraction of atrial myocardium (P wave of ECG). The signal then travels to the AV node where it is delayed, allowing time for the blood to leave the atria before ventricular contraction begins. After the AV nodes, the signal is conducted down the left and right bundle branches and the Purkinje fibers to the ventricular myocardium. Depolarization of the ventricles is recorded as the QRS complex in the ECG, and repolarization of the ventricles is recorded as the T wave. This process is summarized in Fig. 1. In the last lab, we saw how the heart's electrical activity can be analyzed using an ECG (specifically, the standard bipolar lead II). This week, we will explore how to interpret the magnitude and direction of the ECG signal and how to estimate the mean electrical axis of the QRS complex.

Einthoven's Law

ECGs are produced by a system of leads attached to the body. One common lead system uses the three standard bipolar limb leads. The electrode placement defines the recorded direction of the lead, when going from the negative to the positive electrode. The ECG is read as the computed difference (magnitude) between the positive and negative electrodes as a change in voltage over time.

The standard bipolar limb leads, their polarity, and axes are as follows:

Lead	Polarity	Lead Axis
Lead I	right arm (-) to left arm (+)	180° to 0°
Lead II	right arm (-) to left leg (+)	-120° to 60°
Lead III	left arm (-) to left leg (+)	-60° to 120°

The placement of the leads is the result of the extension of Einthoven's Triangle, a triangle of three leads that surrounds the heart. The ECG for each lead displays the amplitude of the voltage as projected onto that lead. To determine the value closest to the actual direction of the heart, Lead II is usually used. However, because the three leads are connected in a loop, Kirchoff's Voltage Law (KVL) states that the voltage drop around the loop must be equal to zero, and therefore, Einthoven's Law can be used to determine the value of one lead if we know the value of the other two.

Mean Electrical Axis: Direction

During the cardiac cycle, the current spreads along specialized pathways, depolarizing the heart in a specific sequence. Consequently, the electrical activity has directionality, that is, a spatial orientation represented by an electrical axis. The summation of all the vectors occurring in a cardiac cycle (i.e. the preponderant direction of current flow during the cardiac cycle) is called the mean electrical axis (Fig. 2)

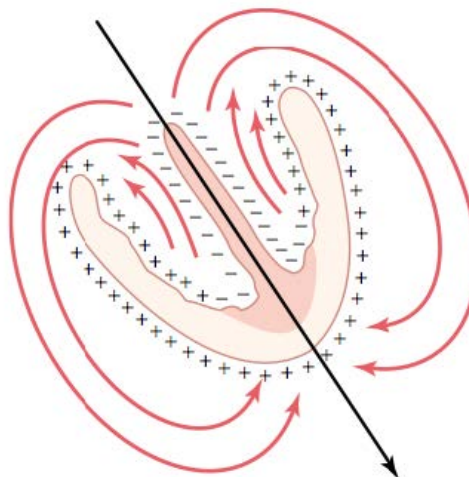


Figure 2: Mean vector through partially depolarized ventricles

Since the QRS complex caused by ventricular depolarization represents the majority of the electrical activity of the heart, you can approximate the mean electrical axis by looking only in this interval. Typically, in an adult, the mean electrical axis lies along a

line extending from the base of the apex of the heart and to the left of the interventricular septum pointing toward the lower left rib cage.

Mean Electrical Axis: Magnitude

The magnitude of the recorded voltage in the ECG is directly proportional to the amount of tissue being depolarized. Most of the mass of the heart is made up of ventricular myocardium. Therefore, the largest recorded waveform, the QRS complex, reflects the depolarization of the ventricles. Furthermore, since left ventricular mass is significantly greater than right ventricular mass, more of the QRS complex reflects the depolarization of the left ventricle and thus the orientation of the mean electrical axis is to the left of the interventricular septum.

ECG Interpretation

A good mathematical tool for representing the measurement of a lead is the vector. At any given moment during the cardiac cycle, a vector may represent the net electrical activity seen by a lead (Fig. 3). The length of the arrow represents the magnitude of the electrical current. The orientation of the arrow represents the direction of current flow. The tip of the arrow represents the positive pole of the electrical current, and the tail represents the negative pole.

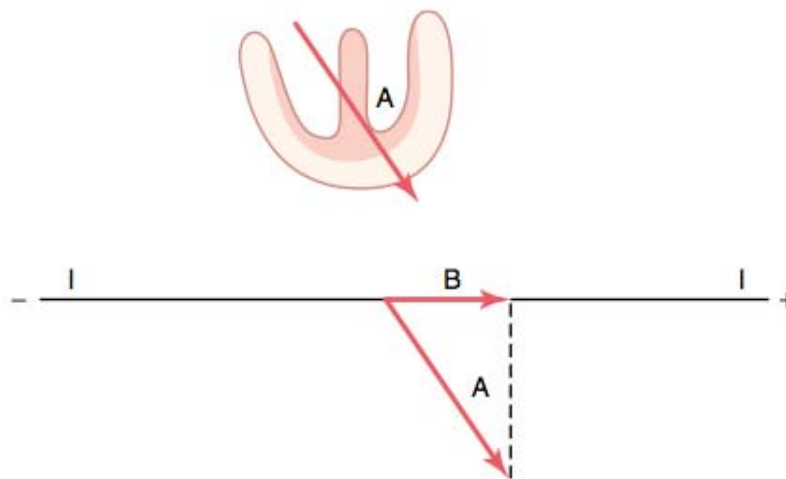


Figure 3: Projection of the mean electrical axis onto lead I (vector B)

Experimental Methods

Electrode Placement

1. Place electrodes following the configuration for lead I:
 - Red (+): left arm
 - White (-): right arm
 - Black (ground): right ankle

2. Place electrodes following the configuration for lead III:
 - Red (+): left leg
 - White (-): left arm
 - Black (ground): right ankle
3. Plug leads I and III into channels 1 and 3, respectively.
4. Open BIOPAC student lab lessons software and select *Create/Record a new experiment*, and underneath it, *Create empty graph* (Fig. 4).

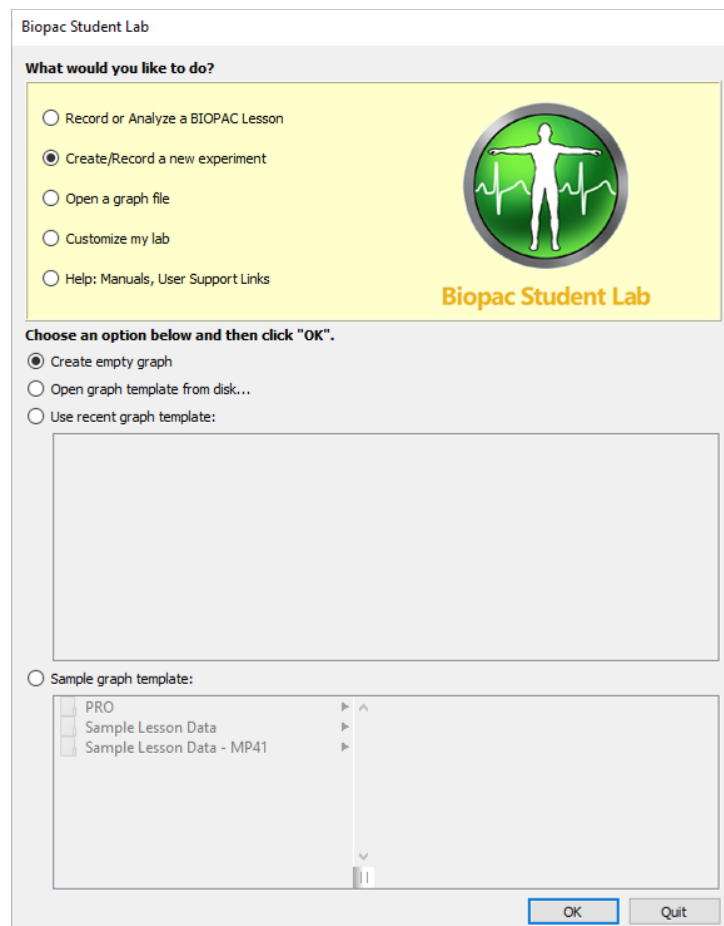


Figure 4: Creating a new BIOPAC experiment

5. Select *Electrocardiogram (ECG)*, 0.05 - 150 Hz for channels 1 and 3 under presets, then check all the boxes for channels 1 and 3 (Fig. 5). Click *Close* when you are done.

Acquire	Plot	Value	Channel	Label	Preset	Channel Sampling Rate
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	CH1	Electrocardiogram (ECG)	Electrocardiogram (ECG), .05	250.000 Hz
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	CH2	Analog input	none	2.000 kHz
<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	CH3	Electrocardiogram (ECG)	Electrocardiogram (ECG), .05	250.000 Hz
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	CH4	Analog input	none	2.000 kHz

Figure 5: Select presets and check all the boxes for channels 1 and 3

Data Collection

You will perform data acquisition for four different conditions (supine, sitting up, inhaling, and exhaling). Calibration is not needed for this lab. Remember to insert markers to separate conditions, and to save the data file at the end.

Supine

1. Have the subject lay in supine position on a mat with eyes closed. Make sure the subject remains relaxed and still.
2. Collect data for 20 seconds (press the green *Start* button on the upper left hand corner).
3. Click the *Stop* button to pause or stop recording.

Sitting Up

1. Have the subject sit upright in a chair. Make sure the subject remains relaxed and still.
2. Collect data for 10 seconds.
3. Click the *Stop* button to pause or stop recording.

Inspiration and Expiration

1. Have the subject remain seated.
2. Collect data while the subject deeply inspires and expires, approximately 5 times in a period of 20 seconds.
3. Place a marker after each inspiration and expiration, making sure to label them properly.
4. Click the *Stop* button to pause or stop recording.

Data Analysis

1. Record the maximum electrical magnitude (amplitude) of the R waves for each condition for each of the two leads in the handout. Make sure the data corresponds to the correct lead number.
2. Using data from the table, graphically determine the Mean Electrical Magnitude and Mean Electrical Axis for the conditions of Supine and Sitting Up. You may plot the values from both conditions on the graph on the handout.
3. Explain the difference (if any) in the amplitudes of Leads I and III, as well as the Mean Electrical Magnitude and Axis under the two conditions.
4. Using data from the table, graphically determine the Mean Electrical Magnitude and Mean Electrical Axis for the conditions of Inhaling and Exhaling. You may plot the values from both conditions on the graph on the handout.

5. Explain the difference (if any) in the amplitudes of Leads I and III, as well as the Mean Electrical Magnitude and Axis under the two conditions.
6. Give normal ranges of the mean electrical axis. Is your data within range?
7. What factors affect the orientation of the mean electrical axis (list at least 2)?
8. What factors affect the amplitude of the R wave recorded on the different leads (list at least 2)?
9. It is unlikely your ECG data looks exactly like "textbook" data. List at least 3 sources of experimental error.